# GESKE

GERMAN BEAUTY TECH

## Quickstart Guide

[www.geske.com](http://www.geske.com)

SmartAppGuided™  
**Lip Volumizer & Booster | 6 in 1**


GESKE Beauty Tech GmbH,  
Leipziger Platz 18, 10117 Berlin, Germany  
[info@geske.com](mailto:info@geske.com), [www.geske.com](http://www.geske.com)

# Short Introduction



Scan for the full instructions  
in your language

Soft Pulse Instant  
Lip Volumizing

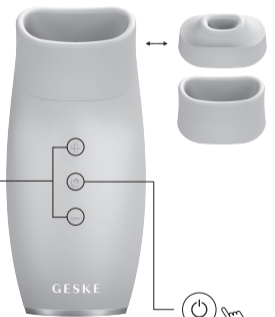Entire Lip Cap  
On Spot Lip Cap

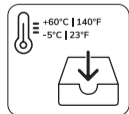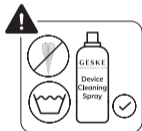GESKE German Beauty Tech App  
Available in 45 Languages

Apple and the Apple logo are trademarks of Apple Inc.  
Google Play and the Google Play logo are trademarks of Google LLC.

**The safety information must also be read before first use. Keep both instructions for use.**

## YOUR PRODUCT AT A GLANCE

Your product is a multi-functional lip plumper. It is safe to be used up to 2 times a day, please do not use more frequently. For maximal safety and ideal results, always follow your personalized Skin Care Routine created by the GESKE App.

**Download the GESKE App and get your personalized routine.**

### SPECIFICATIONS:

**Material:** ABS, Silicone | **Storage temperature:** -5 - +60 °C/23-140 °F | **Battery:** Li-Ion 3.7V | **Battery Capacity:** 380 mAh | **Maximum Charging Current:** 320 mA | **Charging Voltage:** 5V | **Charging Time:** approx. 2,5 hours | **Operating Time:** Operating time differs depending on activated technologies & intensities | **Interface:** USB | **Programs:** 3 Suction Modes with different suction and pause times

### IMPORTANT SAFETY INFORMATION! READ BEFORE USE!

**The product is not suitable for use under the following conditions:**

- Start with the lowest intensity and slowly increase as needed
- Never use the suction function in the same place longer than recommended by the GESKE app to avoid irritation or slight bruising
- Do not use the product on irritated or injured lips or open wounds (after any kind of injury), in the case of bacterial or viral infections of the skin (either in the lip area and elsewhere)
- Do not use the product if it is 1) untypically warm, 2) mechanically damaged or deformed, 3) discolored
- If pain and/or discomfort occurs during use, discontinue use and seek medical advice
- Only use the product as recommended in your personalized routine
- For hygiene reasons, we recommend you do not share your product with anyone else
- The product is intended for cosmetic use on the lips, GESKE assumes no liability for harmful consequences arising from the use on other body areas or other improper use
- Always make sure the plug and socket are completely dry before charging. Otherwise, electric shocks, a short circuit or fire may result. The product must not be used while connected to the power outlet

**The product is not suitable for use under the following conditions:**

- If you suffer from hemorrhagic diseases (congenital blood clotting disorders) and/or take medication that affects blood clotting
- If you suffer from abnormal changes in the facial areas or herpes labialis

- Shortly and in the first weeks after cosmetic lip surgery or lip injections
- With tangible impurities or unnaturally hardened areas in the lip

**In the following cases, a consultation with a physician is necessary as to whether the product is suitable for you to use:**

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